



**NEW**  
**OWNER'S MANUAL,**  
**TUBE INSTALLATION**  
**AND**  
**OPERATING**  
**INSTRUCTIONS**  
**FOR ALL**  
**SPERTI SUNLAMPS**

*FOR PROPER USAGE, PLEASE READ  
YOUR OWNER'S MANUAL THOROUGHLY BEFORE  
STARTING YOUR SUNLAMP.  
THE LABEL ON YOUR SUNLAMP ALSO CONTAINS  
PROPER USAGE INSTRUCTIONS.*

Keep the carton for storage and as a shipping container  
in the event service is required. Original packaging  
is required for all warranty work.

**TOLL FREE #**  
**1-800-544-3757**

## TABLE OF CONTENTS

	<b>Page</b>
Safety Precautions .....	1
Safety Precautions/Tanning .....	2
Skin Type / Recommended Exposure .....	3
Operating Instructions (All Models) .....	4
Operating Instructions (contd) Caution/Improper Handling	5
Trouble Shooting .....	6
Trouble Shooting Guide .....	7
Ultraviolet Tube Facts .....	8
Back Labels .....	9
Warranty & Service Information .....	10

Please take a moment and examine your Sunlamp. Make note below of the Model Number, Wattage, and Date of Purchase. Keep this Manual handy when contacting us for service or further information.

**PROOF OF PURCHASE REQUIRED FOR  
ALL WARRANTY WORK.**

You can be sure of the quality of your **Sperti Suntube** - we have been making them since **1940**. Our sunlamps are in conformity with performance standards for sunlamp products under 21 CFR part 1040.20. Your sunlamp has an automatic shut-off timer. Your **SPERTI SUNLAMP** can do everything the sun can do, and more. It can tan you at night, tan you during the day, tan you indoors, and turn itself off.

Included with your Sperti Sunlamp: Protective eyewear (1), instructions, warranty/registration cards

### *FOR FURTHER INFORMATION OR QUESTIONS*

*WRITE:* KBD, Inc.  
2550 American Court  
Crescent Springs, KY 41017 USA  
859-331-0800  
Fax 859-331-0802  
sales@sperti.com  
www.sperti.com

## **SAFETY PRECAUTIONS:**

- DO consult a physician before using a sunlamp if you have a skin disease, are taking medication or pregnant
- DO wear protective eyewear when using a sunlamp to avoid injury to the eyes. Protective eyewear must be worn by all persons in a room in which a sunlamp is in operation. Additional eyewear may be purchased by contacting your distributor or **SPERTI SUNLAMP**, 2550 American Court, Crescent Springs, KY 41017, 859-331-0800, sales@sperti.com
- DO unplug sunlamp before cleaning or servicing your sunlamp.
- DO follow recommended distance and exposure times. Refer to chart on Page 3.
- DO allow your sunlamp to cool to room temperature (10 to 15 minutes) before restarting sunlamp. This cooling period must be allowed even when your sunlamp fails to advance into the blue glow.
- DO operate sunlamp at room temperature. Dampness, or humidity, may have some effect on the operation of your lamp. Blue glow will appear, and then bright white after reaching full power. Reset timer to desired exposure time.
- DO call, write, fax or E-mail Sperti Customer Service, 2550 American Court, Crescent Springs, KY 41017, 859-331-0800, 859-331-0802 (fax), sales@sperti.com, if you believe your sunlamp is defective. You should not return your sunlamp to the retail store where it was purchased. Problems are often solved quickly by telephone.

\* \* \* \* \*

- DO NOT fall asleep under a sunlamp.
- DO NOT allow children to operate or handle any sunlamp.
- DO NOT use a sunlamp in a bathroom or near water to avoid the possibility of electrical shock.
- DO NOT leave sunlamp while sunlamp is in operation.
- DO NOT use a sunlamp after it has been damaged or dropped. (Consult our Customer Service Department 859-331-0800, 859-331-0802 (fax), E-mail sales@sperti.com.

- DO NOT use sunlamp as a heating or drying device.
- DO NOT use sunlamp on an infant or invalid.
- DO NOT use sunlamp near combustible materials at any time.
- DO NOT touch metal lamp housing during or immediately after operation of lamp. The metal surface becomes HOT.
- DO NOT block vents on sunlamp.
- DO NOT position sunlamp in a heavily traveled area.
- DO NOT operate sunlamp without wire guard.
- DO NOT MOUNT OR HANG YOUR SUNLAMP ON A WALL, CEILING, DOOR, ETC.
- DO NOT apply lotion or cosmetics to skin. Application of creams, lotion and other cosmetics could have an effect on the skin's reaction to your **SPERTI SUNLAMP**.
- DO NOT point or position a sunlamp directly downward or on its side. Heat build-up can shorten the life of your tube.

## **ABUSE OR MISUSE WILL VOID WARRANTY**

### **TANNING**

*If you cannot tan from the sun, you are unlikely to tan using indoor tanning equipment.*

**Pregnant women should consult with their physician prior to tanning.**

Tanning results from the production of Melanin. Melanin is produced when the skin is exposed over a period of time to ultraviolet light, either from the sun or a **SPERTI SUNLAMP**.

One's ability to develop a tan depends usually on body chemistry. Those who produce a high level of Melanin normally develop a deeper tan than those who produce low levels of Melanin. Therefore, the

rate at which you develop a tan and the degree of your tan will depend a great deal on your skin type, not on the light source.

It is essential you take time to develop a tan with your **Sperti Sun-lamp**. It may be 7 to 10 days before a change is noticed, depending on skin type. Regular use according to the instructions is the key to developing and maintaining a tan.

## RECOMMENDED EXPOSURE SCHEDULE

Different skin types require different exposure times and distances. The figures given for time and distances are recommended for average skin, Type III. Not to exceed 3 exposures per week, during first three weeks.

Initial exposure should not exceed 2 minutes for Type II skin or 3 minutes for Type III, IV or V Skin at a distance of 20 inches per section. Exposure at distance less than 20 inches is not recommended. Distance is determined from metal guard **not** light bulb.

	SKIN TYPE	WEEK 1 3 EXPOSURES	WEEK 2 3 EXPOSURES	WEEK 3 3 EXPOSURES	MAINTENANCE 2 TIMES PER WEEK (MAX)
Sensitive skin, very fair, red hair, dense freckles	I	DO NOT USE THIS APPLIANCE.			
Burns easily, Tans minimally or lightly.	II	2 minutes	4 minutes	6 minutes	10 minutes
Often referred to as "average" complexion. Burns moderately, tans about average.	III	3 minutes	5 minutes	7 minutes	10 minutes
Burns minimally, tans easily and above average with each exposure.	IV	3 minutes	5 minutes	7 minutes	10 minutes
Rarely burns, tans easily and substantially.	V	3 minutes	5 minutes	7 minutes	10 minutes

**NOTE:** A section is approximately 12 inches x 12 inches.

**NOTE:** Reddening is usually not visible until 4-8 hours after exposure.

## OPERATING INSTRUCTIONS (ALL MODELS)

Insert plug in 110-120v, 60/50Hz, AC outlet.



### PLUG SAFETY:

Your sunlamp has a polarized plug (one blade is wider than the other), and will fit your outlet in only one way. If the plug does not fit, reverse the plug. If it still does not fit properly, contact a qualified electrician. NEVER FORCE PLUG OR BLADES INTO OUTLET. DO NOT USE AN EXTENSION CORD WITH YOUR SUNLAMP.

**PROTECTIVE EYEWEAR MUST BE WORN AT ALL TIMES DURING OPERATION OF ALL MODELS.**

### I. Models P181, S363F, P164R AND MARK IV-F

**Set timer to 10 minutes.**

DEPRESS Starter Button and **hold only until coils inside tube turn orange** (usually 4-6 seconds), **then immediately release button**. A blue glow will appear, changing to bright white after lamp reaches full power (approximately 3-4 minutes).

### IMPORTANT

Immediately release the switch or button after coils turn orange. Otherwise, the tube may overheat, and the Sunlamp must then cool for about 10-15 minutes before it can be restarted.

**NOTE:** Refer to Trouble-Shooting Guide if the coils do not turn orange or the blue glow does not appear.

### HEATLAMP MODELS: (HT-10 AND HT-10T, AND MARK IV-F)

Place unit at a distance of 24 inches or more. Measure from wire guard. Depress Heat Switch. If more than 10 minutes is desired, reset timer.

## **II. Models DS2000 (Del Sol) and SS2000 (Sentosa Sun)**

Press **ON/MIN** button once. Repeat this process to advance time to 10 minutes (Maximum). Wait 5 seconds, click indicates the start process, causing the bulb coils to glow orange. In 7 seconds another click will ignite your Sperti Sun Tube. AFTER bulb reaches full power (3-4 minutes) use **ON/MIN** button to reset exposure time.

**OFF** button terminates all functions.

## **III. Floor Models S-363-F, Mark IV-F, Sentosa Sun**

Assemble and attach stand; instructions included with your Sunlamp.

Adjust Sunlamp to desired position. NEVER POINT A STAND MODEL STRAIGHT DOWN OR TURN A TABLE MODEL ON ITS SIDE. Heat buildup can shorten the life of your tube.

**CAUTION: IMPROPER HANDLING (E.G., FAILURE TO PROPERLY PLACE TUBE INSIDE OF CONTACTS) MAY RESULT IN A BROKEN GLASS EXHAUST. THE BREAKAGE THAT OCCURS DURING INSTALLATION IS NOT A WARRANTY ITEM, AND A REPLACEMENT TUBE WILL BE AT CUSTOMER EXPENSE.**

## TROUBLE SHOOTING/LAMP INSTRUCTIONS

(When lamp fails to start)

1. Make sure lamp is disconnected before checking tube.
2. Remove protective wire guard.
3. Remove two screws which secure the tube shields and remove both shields.
4. Sometimes in shipping the tube holders (or terminals) are jarred and not making the proper connection. (See diagram #1)
5. With your thumbs move terminals slightly away from the glass ends of tube to make sure there is a good connection. (See diagram #2)
6. Reinstall tube shields and protective guard.
7. Plug in your lamp and try again.

DIAGRAM # 1




TUBE HOLDERS  
(Terminals)

DIAGRAM # 2

TUBE HOLDERS  
(Terminals)

← Move Terminals  
Slightly  
This Direction

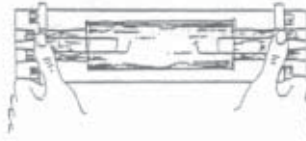

Move Terminals  
Slightly  
This Direction →

**NOTE:** All Sperti Sunlamps automatically shut off. The timer for all models can be seen through the protective eyewear. Once off, the Sunlamp cannot be restarted until it has cooled for approximately 10 to 15 minutes.

## TROUBLE SHOOTING GUIDE

PROBLEM	CORRECTION
<p><b>Coils will not turn orange</b></p>	<ol style="list-style-type: none"> <li>1. No current</li> <li>2. Failure to set timer</li> <li>3. UV Tube not making proper connection</li> </ol>
<p><b>Blue/White glow will not appear</b></p>	<ol style="list-style-type: none"> <li>1. Starter/Sun Button has been held too long and UV Tube has overheated</li> <li>2. UV Tube cracked or faulty</li> <li>3. Sunlamp has not been allowed to cool</li> </ol>
<p><b>Lamp does not tan</b></p>	<ol style="list-style-type: none"> <li>1. Dirty or dusty UV Tube</li> <li>2. Distance from Sunlamp is too great</li> <li>3. Facial cosmetics are blocking UV</li> </ol>
<p><b>Moisture in UV Tube</b></p>	<ol style="list-style-type: none"> <li>1. Mercury Droplets</li> </ol>
<p><b>Black spots on UV Tube</b></p>	<ol style="list-style-type: none"> <li>1. Mercury attached to UV Tube glass</li> </ol>
<p><b>Sunlamp Is smoking and/or giving off odor. (Temporary condition)</b></p>	<ol style="list-style-type: none"> <li>1. Oil on reflector</li> <li>2. Incorrect positioning of Sunlamp</li> </ol>

ALWAYS RETURN YOUR SUNLAMP TO OUR FACTORY FOR REPAIR OR REPLACEMENT, NOT TO THE RETAIL OUTLET WHERE YOUR SUNLAMP WAS PURCHASED.

## ULTRAVIOLET TUBE FACTS



**IMPORTANT NOTE:** There is no connecting filament. Your tube should look like the illustration at left.

- A. The UV Tube has been carefully manufactured to provide many hours of safe tanning. Under normal use, you can expect a life of 250 hours. (10 minutes x 3 sessions per week x 52 weeks = 3 years)
- B. Your UV Tube will last longer if you do not point your Sunlamp directly downward or position it on its side; leave Sunlamp on an angle. Also, wipe your UV Tube with a **dry**, soft towel or cloth **before** plugging in your Sunlamp to receive the maximum benefits from your Tube.

**IMPORTANT.-** If the Sunlamp is allowed to shut off automatically, it must be allowed to cool 10 to 15 minutes before it can be restarted.

**NOTE:** If lamp still fails to start, please contact our Customer Service Department at (859) 331-0800, 800-544-3757, 7:30 AM and 4:00 PM, EST. E-mail [sales@sperti.com](mailto:sales@sperti.com).

---

### OPERATING INSTRUCTIONS FOR YOUR SPERTI SUNLAMP

- Plug cord into 110-120V, 50/60 Hz A.C. receptacle.
- Be sure to put on your protective eyewear enclosed with every lamp.
- Press On/Min button, release. Using top button, select SUN. Set timer to 10 mins., lamp will automatically ignite in 12 seconds.
- A "Warm Up" period of 3-5 minutes is required before you begin your tanning session. After the warm-up period, the timer can be re-set to desired exposure.
- The "OFF" button terminates ALL functions. Your lamp will not re-ignite for 10 minutes.
- RADIANT/HEAT Press On/Min button. Using top button, select HEAT. Set timer.
- After exposure unplug the lamp.

**MANUFACTURED S**



Model Nos. SUNCRUZ  
No. (SC2002B) Floor Stand  
(SC2002B-T) Table Model  
120 VOLT 50/60 HZ. AC ONLY  
800 WATTS.  
Questions? 1-800-544-3757  
Website: [www.sperti.com](http://www.sperti.com)

**RECOMMENDED EXPOSURE** - Initial exposure should not exceed 3 minutes at a distance of not less than 14 inches - 0.36 meter. Use a marked tape to measure the distance. **CAUTION** - Exposure at distances less than 14 inches - 0.36 meter - is not recommended. Different skins require different times and distances of exposure. The figures given above for time and distance are safe for average skin. After recommended exposure, reddening or sunburn is not visible until 4 to 8 hours. Each day increase the exposure time not more than one minute until you have attained the desired degree of tanning. Recommended maximum exposure time is 15 minutes. For maintenance of this degree of tanning, do not exceed your last exposure. Let your skin guide you. Exposure should not be repeated more than once daily. Do not fall asleep under the sunlamp.

**RADIANT HEAT TREATMENT:** Set time to desired time and switch to HEAT.

**TO REPLACE SPERTI SUNLAMP ULTRAVIOLET TUBE, UV - 93** - The tube ends are covered with metal shields. WITH LAMP UNPLUGGED, remove wire guide and the two screws which hold these shields. Remove the shields. Remove the tube with a turning motion. Insert the new tube with red mark on your right and facing you. Rotate tube into holder. Replace the shields and screws and wire guard.

**SPERTI SUNLAMP, KBD, INC., - CRESCENT SPRINGS, KENTUCKY 41017**

**CAUTION** - PROTECT YOUR EYES - WEAR PROTECTIVE EYEWEAR - SPERTI CAT. NO. 21050A. PROTECTIVE EYEWEAR MUST BE WORN BY ALL PERSONS IN ROOM WITH AN OPERATING LAMP TO AVOID INJURY TO THE EYES. READ THE INSTRUCTION BOOKLET BEFORE USING LAMP. DO NOT USE IN BATHROOM. This product complies with 21 CFR 1040.20.

**DANGER: ULTRAVIOLET RADIATION. FOLLOW INSTRUCTIONS. AVOID OVEREXPOSURE. AS WITH NATURAL SUNLIGHT, OVEREXPOSURE CAN CAUSE EYE AND SKIN INJURY AND ALLERGIC REACTIONS. REPEATED EXPOSURE MAY CAUSE PREMATURE AGING OF THE SKIN AND SKIN CANCER. WEAR PROTECTIVE EYEWEAR; FAILURE TO MAY RESULT IN SEVERE BURNS OR LONG-TERM INJURY TO THE EYES. MEDICATIONS OR COSMETICS MAY INCREASE YOUR SENSITIVITY TO ULTRAVIOLET RADIATION. CONSULT PHYSICIANS BEFORE USING SUNLAMP IF YOU ARE USING MEDICATIONS OR HAVE A HISTORY OF SKIN PROBLEMS OR BELIEVE YOURSELF ESPECIALLY SENSITIVE TO SUNLIGHT. IF YOU DO NOT TAN IN THE SUN, YOU ARE UNLIKELY TO TAN FROM THE USE OF THIS PRODUCT.**

**OPERATING INSTRUCTIONS FOR YOUR SPERTI SUNLAMP**

- Wear protective eyewear (provided with unit)
- Examine distance 20 inches
- Adjust lamp position. Avoid pointing directly down. Head bulb on shorter tube lit.
- Plug cord into 110-120V, 60/50 Hz grounded receptacle only.
- To start, set timer to 10 minutes.
- Push and hold button switch. After 2 to 3 seconds, orange glow will appear, release button switch, tube will light with a pale blue color.
- 2 minutes, the tube will glow bright white. Reset timer for desired exposure.
- Unplug lamp for storage and service.

**RECOMMENDED EXPOSURE:** Refer to Exposure Label or Owner's Manual.

**TO REPLACE SPERTI SUNLAMP ULTRA VIOLET TUBE UY - 81** - The tube ends are covered with metal shields. WITH LAMP UNPLUGGED, remove the two screws which hold these shields. Remove the shields. Gently spread the metal strips which hold the tube in place. Remove the tube with a turning motion. Again, gently spread the metal strips and insert the new tube with red mark on your right and facing you using a turning motion. Replace the shields and screws.

**CAUTION - PROTECT YOUR EYES - WEAR PROTECTIVE EYEWEAR - SPERTI CAT. No. 21050A. PROTECTIVE EYEWEAR MUST BE WORN BY ALL PERSONS IN ROOM WITH A OPERATING LAMP TO AVOID INJURY TO THE EYES. READ THE INSTRUCTION BOOKLET BEFORE USING LAMP. DO NOT USE IN BATHROOM.** This product complies with 21 CFR 1040.20.

**DANGER: ULTRAVIOLET RADIATION. FOLLOW INSTRUCTIONS. AVOID OVEREXPOSURE. AS WITH NATURAL SUNLIGHT, OVEREXPOSURE CAN CAUSE EYE AND SKIN INJURY AND ALLERGIC REACTIONS. REPEATED EXPOSURE MAY CAUSE PREMATURE AGING OF THE SKIN AND SKIN CANCER. WEAR PROTECTIVE EYEWEAR. FAILURE TO MAY RESULT IN SEVERE BURNS OR LONG-TERM INJURY TO THE EYES. MEDICATIONS OR COSMETICS MAY INCREASE YOUR SENSITIVITY TO ULTRAVIOLET RADIATION. CONSULT PHYSICIAN BEFORE USING SUNLAMP IF YOU ARE USING MEDICATIONS OR HAVE A HISTORY OF SKIN PROBLEMS OR BELIEVE YOURSELF ESPECIALLY SENSITIVE TO SUNLIGHT. IF YOU DO NOT TAN IN THE SUN, YOU UNLIKELY TO TAN FROM THE USE OF THIS PRODUCT.**

KBD, INC., CRESCENT SPRINGS, KY 41017 1-800-544-3757 www.sperti.com



Model No. P-124-B Table Model  
Model No. S-883-P Floor Model  
AC ONLY 120 volts, 60/50 Hz 500 watts

0270  
Rev. 10/08

**OPERATING INSTRUCTIONS FOR YOUR SPERTI SUNLAMP J F M A M J J A S O N D 05 06 07 08 09**

- Wear protective eyewear (provided with unit)
- Examine distance 20 inches
- Refer to owners manual for exposure schedule
- Adjust lamp position. Avoid pointing directly down. Head bulb on shorter tube lit.
- Plug cord into 110-120V, 60/50 Hz grounded receptacle only.
- To start, set timer to 10 minutes.
- Push and hold button switch. After 2 to 3 seconds, orange glow will appear, release button switch, tube will light with a pale blue color.
- 2 minutes, the tube will glow bright white. Reset timer for desired exposure.
- Unplug lamp for storage and service.

**RECOMMENDED EXPOSURE:** Refer to owners manual. Initial exposure should not exceed 2 minutes at a distance of not less than 20 inches. Subsequent not to exceed 3 times weekly. After recommended exposure, recharging or rearming usually is not needed and 4 to 6 hours. Recommended maximum exposure times is 15 minutes. For maintenance of skin, do not exceed your set exposure. Let your skin guide you. Do not feel exposed upon the skinning.

**TO REPLACE SPERTI SUNLAMP ULTRA VIOLET TUBE UY - 81** - The tube ends are covered with metal shields. WITH LAMP UNPLUGGED, remove the two screws which hold these shields. Remove the shields. Gently spread the metal strips which hold the tube in place. Remove the tube with a turning motion. Again, gently spread the metal strips and insert the new tube with red mark on your right and facing you using a turning motion. Replace the shields and screws.

**CAUTION - PROTECT YOUR EYES - WEAR PROTECTIVE EYEWEAR - SPERTI CAT. No. 21050A. PROTECTIVE EYEWEAR MUST BE WORN BY ALL PERSONS IN ROOM WITH A OPERATING LAMP TO AVOID INJURY TO THE EYES. READ THE INSTRUCTION BOOKLET BEFORE USING LAMP. DO NOT USE IN BATHROOM.** This product complies with 21 CFR 1040.20.

**DANGER: ULTRAVIOLET RADIATION. FOLLOW INSTRUCTIONS. AVOID OVEREXPOSURE. AS WITH NATURAL SUNLIGHT, OVEREXPOSURE CAN CAUSE EYE AND SKIN INJURY AND ALLERGIC REACTIONS. REPEATED EXPOSURE MAY CAUSE PREMATURE AGING OF THE SKIN AND SKIN CANCER. WEAR PROTECTIVE EYEWEAR. FAILURE TO MAY RESULT IN SEVERE BURNS OR LONG-TERM INJURY TO THE EYES. MEDICATIONS OR COSMETICS MAY INCREASE YOUR SENSITIVITY TO ULTRAVIOLET RADIATION. CONSULT PHYSICIAN BEFORE USING SUNLAMP IF YOU ARE USING MEDICATIONS OR HAVE A HISTORY OF SKIN PROBLEMS OR BELIEVE YOURSELF ESPECIALLY SENSITIVE TO SUNLIGHT. IF YOU DO NOT TAN IN THE SUN, YOU UNLIKELY TO TAN FROM THE USE OF THIS PRODUCT.**

KBD, INC., CRESCENT SPRINGS, KY 41017 1-800-544-3757 www.sperti.com



Model P121  
AC ONLY 120 volts, 60/50 Hz 700 watts

0271  
Rev. 10/08

**OPERATING INSTRUCTIONS FOR YOUR SPERTI SUNLAMP**

- Wear protective eyewear (provided with unit)
- Examine distance 20 inches
- Adjust lamp position. Avoid pointing directly down. Head bulb on shorter tube lit.
- Plug cord into 110-120V, 60/50 Hz grounded receptacle only.
- To start, set timer to 10 minutes.
- Push and hold button switch. After 2 to 3 seconds, orange glow will appear, release button switch, tube will light with a pale blue color.
- 2 minutes, the tube will glow bright white. Reset timer for desired exposure.
- Unplug lamp for storage and service.

**RECOMMENDED EXPOSURE:** Refer to Exposure Label or Owner's Manual.

**TO REPLACE SPERTI SUNLAMP ULTRA VIOLET TUBE UY - 81** - The tube ends are covered with metal shields. WITH LAMP UNPLUGGED, remove the two screws which hold these shields. Remove the shields. Gently spread the metal strips which hold the tube in place. Remove the tube with a turning motion. Again, gently spread the metal strips and insert the new tube with red mark on your right and facing you using a turning motion. Replace the shields and screws.

**CAUTION - PROTECT YOUR EYES - WEAR PROTECTIVE EYEWEAR - SPERTI CAT. No. 21050A. PROTECTIVE EYEWEAR MUST BE WORN BY ALL PERSONS IN ROOM WITH A OPERATING LAMP TO AVOID INJURY TO THE EYES. READ THE INSTRUCTION BOOKLET BEFORE USING LAMP. DO NOT USE IN BATHROOM.** This product complies with 21 CFR 1040.20.

**DANGER: ULTRAVIOLET RADIATION. FOLLOW INSTRUCTIONS. AVOID OVEREXPOSURE. AS WITH NATURAL SUNLIGHT, OVEREXPOSURE CAN CAUSE EYE AND SKIN INJURY AND ALLERGIC REACTIONS. REPEATED EXPOSURE MAY CAUSE PREMATURE AGING OF THE SKIN AND SKIN CANCER. WEAR PROTECTIVE EYEWEAR. FAILURE TO MAY RESULT IN SEVERE BURNS OR LONG-TERM INJURY TO THE EYES. MEDICATIONS OR COSMETICS MAY INCREASE YOUR SENSITIVITY TO ULTRAVIOLET RADIATION. CONSULT PHYSICIAN BEFORE USING SUNLAMP IF YOU ARE USING MEDICATIONS OR HAVE A HISTORY OF SKIN PROBLEMS OR BELIEVE YOURSELF ESPECIALLY SENSITIVE TO SUNLIGHT. IF YOU DO NOT TAN IN THE SUN, YOU UNLIKELY TO TAN FROM THE USE OF THIS PRODUCT.**

KBD, INC., CRESCENT SPRINGS, KY 41017 1-800-544-3757 www.sperti.com



Model No. DE, BCL  
Model No. S8170GA SLN  
AC ONLY 120 volts 60/50 Hz 800 watts

0272  
Rev. 05/05

**J F M A M J J A S O N D**  
**RECOMMENDED EXPOSURE**

Different skin types require different exposure times and distances. The figures given for time and distances are recommended for average skin, Type II. Not to exceed 3 exposures per week.

Initial exposure should not exceed 2 minutes for Type II skin or 3 minutes for Type III, IV or V skin at a distance of 20 inches per session. Exposure at distance less than 20" is not recommended. Distance is determined from metal guard set light bulb.

**NOTE:** A session is approximately 12 inches by 12 inches.

**NOTE:** Redding or sunburn is usually not visible until 4-8 hours after exposure.

05 06 07 08 09



	SKIN TYPE	WEEK 1 3 EXPOSURES	WEEK 2 3 EXPOSURES	WEEK 3 3 EXPOSURES	MAINTENANCE 2 TIMES PER WEEK MAX	
Sensitive skin, very fair, red hair, dense freckles	I	DO NOT USE THIS APPLIANCE.				
Burns easily, Tans minimally or lightly.	II	2 minutes	4 minutes	6 minutes	10 minutes	
Often referred to as "average" complexion. Burns moderately, tans about average.	III	3 minutes	5 minutes	7 minutes	10 minutes	
Burns minimally, tans easily and above average with each exposure.	IV	3 minutes	5 minutes	7 minutes	10 minutes	
Rarely burns, tans easily and substantially.	V	3 minutes	5 minutes	7 minutes	10 minutes	

0273  
Rev. 04/02

## LIMITED WARRANTY

### Sperti Sunlamp

**What Is Covered:** All parts defective in materials or quality of workmanship.

**For How Long:** Three (3) Months from date of purchase.

**What We Will Do:** We will decide either to replace or repair free-of-charge any defective part.

**What We Will Not Do:** Pay shipping or transportation costs.

**What You Must Do:** The warranty registration card must be completed within ten (10) days from the date of purchase. The completion and return of the warranty registration card to Sperti Sunlamp, 2550 American Court, Crescent Springs, KY 41017, within ten (10) days of purchase of this Sperti Sunlamp is a condition of warranty coverage.

**How To Get Warranty Services:** For service, defective sunlamp or defective part of the sunlamp must be returned to the Factory Service Department, Sperti Sunlamp, 2550 American Court, Crescent Springs, KY 41017.

**Other Conditions:** This warranty covers the first purchaser only. This warranty does not cover commercial use. This warranty shall only apply to Sperti Sunlamps that have not been tampered with or altered in any way or subjected to misuse, neglect or accident.

**State Law Rights:** Implied warranties are limited to the duration of three (3) months.

Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.

This warranty gives you specific legal rights and you may also have other rights which vary from state to state.



**ALL SPERTI SUNLAMPS  
ARE UL LISTED**